


Enteral Nutrition Administration Omni Pump

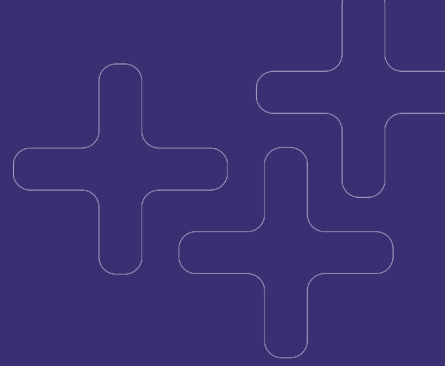

*This information is not intended to be a substitute for professional advice.
As such, please follow recommendations from your doctor, dietitian, nurse or pharmacist.*

Getting started

Use a clean work area. Clean your hands thoroughly using soap and warm water for at least 20 seconds. If no soap and water is available, use an alcohol-based hand sanitizer. Dry your hands with a clean towel.

1. Prepare home tube feeding

- A. Gather your formula and supplies, placing on the clean work area.
(Fig. 1)
- B. Check expiration date on formula carton and make sure the formula name matches your order. Wipe top with a clean cloth and shake well before opening.



Opened cartons of formula should be immediately covered and refrigerated. Product should not hang in the feeding bag for more than 12 hours (if a reconstituted powder is used, it should not hang in the bag for more than 4 hours). Use or discard any opened formula within 24 hours or per manufacturer recommendations.

- C. Place a towel on your lap as formula may spill during administration.
- D. Sit in an upright position or elevate the head of your bed at least 30-45 degrees. **DO NOT** lay flat during feeding and for at least one hour after.

Your Option Care Health team is available
on call 24/7/365 to assist you.


2. Administer tube feeding

- A. Holding the feeding bag upright, open the lid and fill with ordered amount of formula, securely closing the cap afterward. (Fig. 2)
- B. Turn the pump on by pressing the power key on the lower part of the front of the pump. (Fig. 3)
- C. Press the corresponding menu selection button ► **'KEEP SETTINGS'** or ► **'CLEAR SETTINGS'** to perform the desired action.

*Select **'KEEP SETTINGS'** for day to day operation.*

- D. Load the feeding cassette onto the pump by following the animated instructions on the screen. You will hear a “click” when the cassette has been installed correctly. (Fig. 4) The Set Loaded screen will confirm that the feeding cassette has been loaded correctly. Remove the clear plastic cover from the end of the tubing. (Fig. 5)
- E. When ready to prime, press ► **'PRIME'** and then ► **'AUTO PRIME'** to start priming. (Fig. 6)
- F. Press and hold ► **'PRIME FEED'** to manually run the formula to the end of the tubing and then press ► **'DONE.'**
- G. Press ► **'ADJUST RATE'** and then use the + and - buttons to increase or decrease the rate. Select ► **'OK'** when the rate is set. (Fig. 6)
- H. If your feeding tube has a clamp, make sure it is in the closed position before opening the end of your feeding tube. If your tube does not have a clamp, pinch tubing before opening the end.
- I. Open your feeding tube cap, unclamp tubing and flush feeding tube slowly with water. (Fig. 7)
- J. Connect the end of the primed feeding bag to your feeding tube. Press ► **'START'** to begin feeding.
- K. Once feeding is complete, disconnect the feeding bag and use the syringe to flush feeding tube slowly with water.
- L. Close the end of the feeding tube and clamp tubing.
- M. Press the power button to turn pump off.
- N. Rinse feeding bag thoroughly with warm water between feedings.

Use one feeding bag every 24 hours.

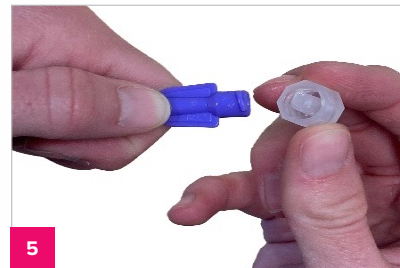

Feeding schedule

Formula name(s):	
Goal amount/day	
Goal rate:	
Interval/frequency:	
Additional notes:	

Administering water flushes

Water flushes are administered to provide adequate hydration and to prevent tube clogging.

- Use a 60 mL syringe to flush your feeding tube slowly with lukewarm water
- Flush with ____mL water before and ____mL water after feeding
- Flush with an additional ____mL water each day

Don't use too much force when flushing your tube.

Tips for success

- **DO NOT** give tube feeding if your stomach is upset or if you are vomiting.
- If you develop an uncomfortable feeling of fullness, slow down the feeding or wait one hour before giving more feeding.
- **DO NOT** put anything other than formula or water into the tube unless specifically instructed to do so by your healthcare provider.

**Your Option Care Health team is available
on call 24/7/365 to assist you.**



Scan the QR code to view video instructions
and additional educational resources or
visit optioncarehealth.com/education

